


# Putting PPE ON 1/1



Contact Gown

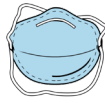

N95  
Respirator



Face shield or  
goggles



Gloves

## Gown

1 **Remove** any personal items and jewelry and put in secure location, not in pockets.

2 **Sanitize** hands.

3 **Put on** contact gown outside room.

Open-end **faces** your back.

**Tie** the back of the gown.

## Mask + Eyes

4 **Put on** gloves.

6 **Fit** mask to nose, make sure no air is escaping.

7 **Put on** face shield or goggles.

8 **Remove** gloves.

5 **Put on** mask.

9 **Sanitize** hands.

## Gloves + Entry

10 **ENTER** room

11 **Sanitize** hands.

12 **Put on** gloves (over the cuffs of the gown if worn).

13 **Do not** touch face or re-adjust mask or face shield inside room.

# taking PPE OFF 1/1



N95, face shield or goggles, gown, gloves

## Gown + Gloves

<p>1 Sanitize gloves.</p>	<p>2 Cross arms and <b>grip</b> gown on shoulders.</p> <p><b>Pull and break</b> gown in controlled fashion.</p>	<p>3 Roll the gown towards your hands. <b>Remove</b> the gloves with the gown.</p> <p><b>Dispose</b> of gloves and gown.</p>	<p>4 Sanitize hands.</p>
---------------------------	---	--	--------------------------

**EXIT patient room**



**EXIT patient room**

## Eyes

<p>5 Sanitize hands.</p>	<p>6 <b>Do not touch face.</b></p>	<p>7 <b>Remove</b> face shield by the <b>strap</b> over your head without touching your skin.</p>
--------------------------	------------------------------------	---

## Mask

<p>8 Sanitize hands again.</p>	<p>9 <b>Pinch</b> bottom strap and <b>pull</b> far over head.</p> <p><b>Do not let</b> straps touch your face.</p>	<p>10 <b>Pinch</b> top strap and <b>pull</b> far over head.</p> <p><b>Do not let</b> straps touch your face as you <b>remove</b> the N95.</p>
--------------------------------	--	---

## Wash

<p>11</p>	<p><b>Head immediately</b> to handwashing station.</p> <p><b>Wash hands</b> with soap and water.</p>
-----------	--